



UmanoMedical

THERAPEUTIC SERIES

airnest orange HYBRID™



**Pressure Injury Prevention · Hybrid System for Treatment Flexibility  
· Pressure Redistribution**

## KEY BENEFITS

- Offers a non-powered pressure redistributing surface with the ability to convert it to a powered surface and offer additional therapies of alternating pressure (AP) and low air loss (LAL).

### Non-Powered Therapy

- SAT (self-adjusting technology) automatically adjusts internal air pressure in reaction to body movements or repositioning.

### Powered Therapies

- On-Demand Low Air Loss
- Alternating Pressure
- Static Mode



**PRODUCT FEATURES**

**AIR FOAM CELLS**


Provides pressure redistribution and patient comfort.

**OPTIONAL CONTROL UNIT**

Intuitive control panel allows the caregiver to quickly set patient weight, select desired therapies, and determine AP cycle times.

**TOP COVER**

Top cover is low friction and low shear force producing, highly vapor permeable, and is available in nylon or polyurethane stretch material.

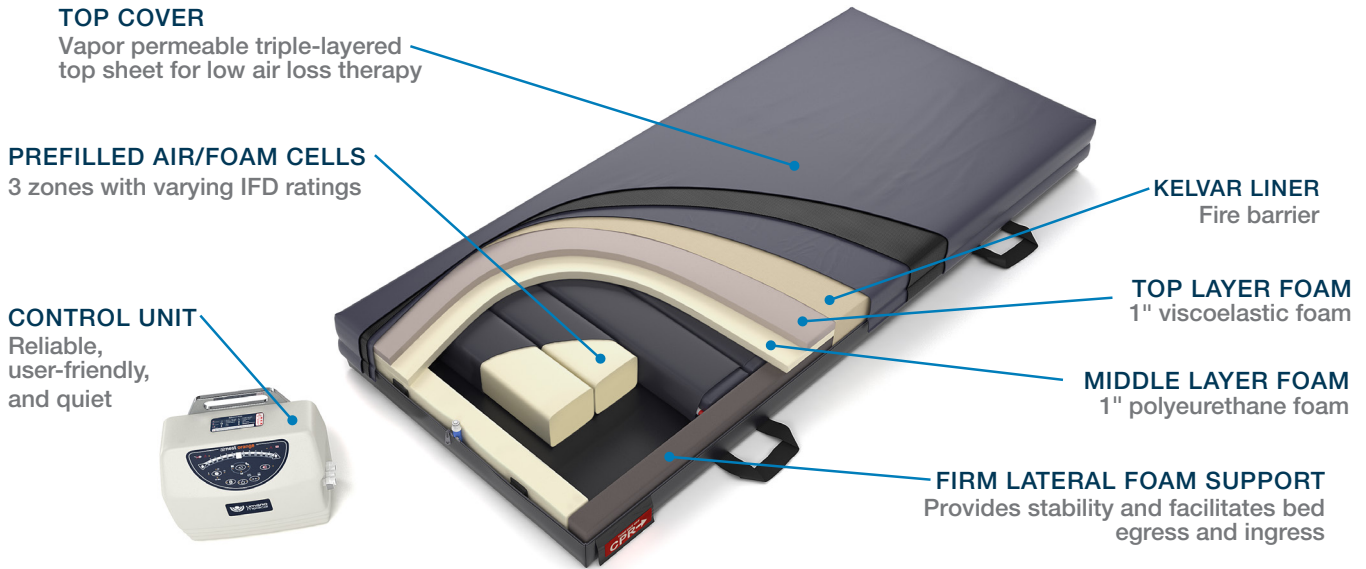

**TECHNICAL SPECIFICATIONS**

**Standard Size** 35" W X 80" or 84" L X 7"H (88.9 cm W X 203 cm or 213 cm L X 17.8 cm H)

**Weight Capacity** 500 lb (227 Kg)

**Mattress Weight** 43.5 lb (19.7 Kg)

<b>Control Unit</b>	Dimensions	13.8" x 11" x 6.5" (35 cm x 28 cm x 16.4 cm)
	Weight	10.2 lb (4.6 Kg)
	Cycle Times	AP Cycle Times: 5, 10, 15, & 20 minutes
	Power Cord	14' hospital grade detachable
	Input Voltage	90 264 V, 50/60 Hz, 1A/0.5A, Universal
	ETL Listed	Conforms to UL STD 60601-1 with respect to Electrical Shock, Fire and Mechanical Hazards
	Certified to	CAN/CSA STD C22.2 No. 601.1

**Flammability** Meets flammability standards CAL 117, CAL 129, 16 CFR 1632, & 16 CFR 1633

Specifications can change without notice.



**Umano Medical Inc.**  
 230, boul. Nilus-Leclerc  
 L'Islet QC G0R 2C0 CANADA  
 Toll Free: 1 (844) 409-4030  
 information@umanomedical.com  
 umanomedical.com

Umano Medical is distributor of the Therapeutic Series manufactured by KAP MEDICAL