



UmanoMedical

COMFORT SERIES

comfortnest blue™




## Pressure Redistribution • Comfort • Durability

### KEY BENEFITS

- 2 layer / 2 zone design offers pressure redistribution and comfort
- Heel zone with low IFD foam for enhanced pressure redistribution in the high-risk heel area
- Entire foam core is constructed of 100% medical grade non-allergenic foam
- All components are anti-microbial and assist in prohibiting growth of bacteria and fungi
- Clam shell design cover features double-folded seams and 4 reinforced handles
- 450 lb weight capacity
- comfortnest blue+ model offers High Resilience (HR) foam and a weight capacity of 950 lb

PRODUCT FEATURES

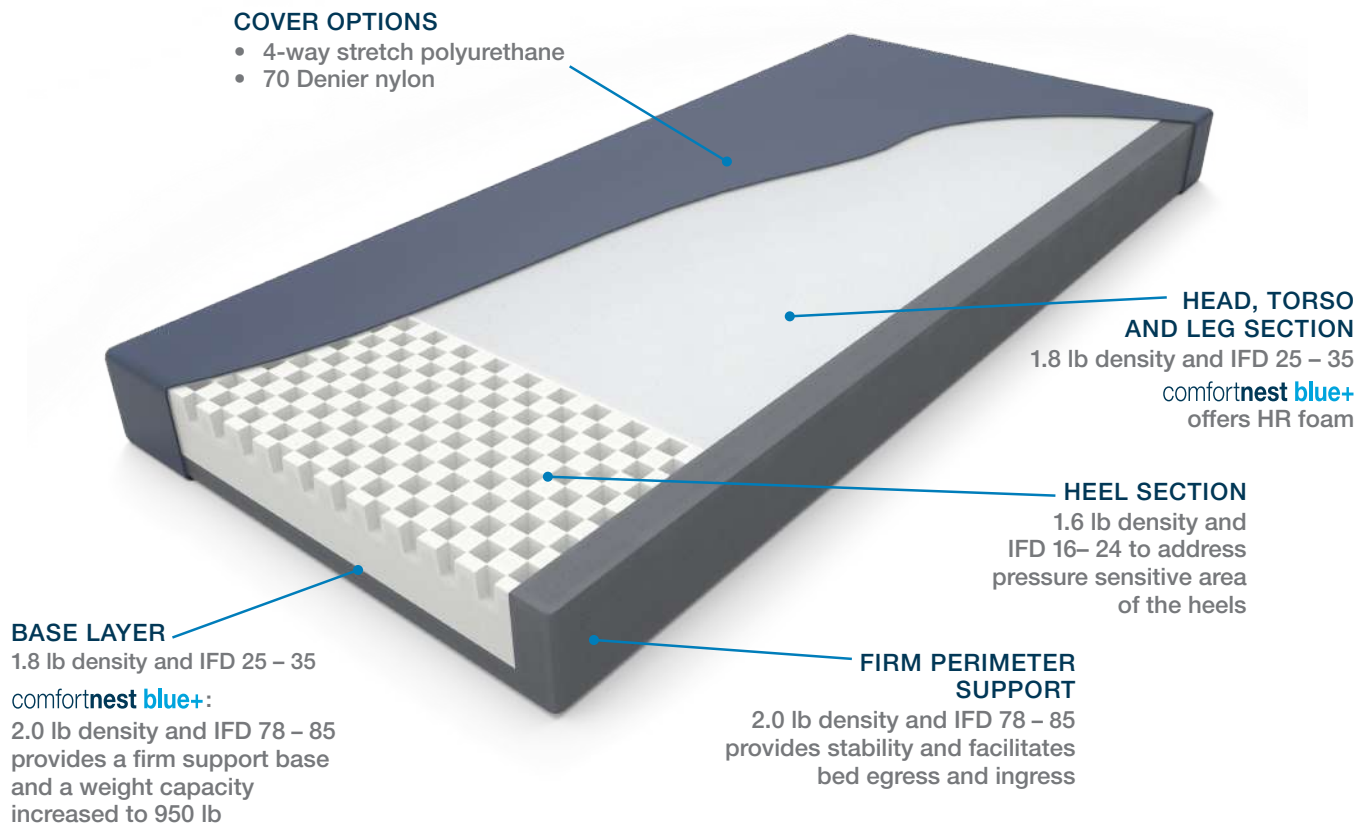

**COVER OPTIONS**

- 4-way stretch polyurethane
- 70 Denier nylon

**HEAD, TORSO AND LEG SECTION**  
1.8 lb density and IFD 25 – 35  
comfortnest blue+ offers HR foam

**HEEL SECTION**  
1.6 lb density and IFD 16– 24 to address pressure sensitive area of the heels

**BASE LAYER**  
1.8 lb density and IFD 25 – 35  
comfortnest blue+:  
2.0 lb density and IFD 78 – 85 provides a firm support base and a weight capacity increased to 950 lb

**FIRM PERIMETER SUPPORT**  
2.0 lb density and IFD 78 – 85 provides stability and facilitates bed egress and ingress

TECHNICAL SPECIFICATIONS

<b>Dimensions</b>	35" W x 80" or 84" L x 6" H (88.9 cm W x 203.2 cm or 213.4 cm L x 15.2 cm H) Custom sizes available
<b>Specialty Model</b>	Behavioral & Mental Health model available: comfortnest blue <sup>MH</sup> Fully enclosed cover, no handles, no flaps, no zippers.
<b>Weight Capacity</b>	450 lb (205 Kg) 950 lb (432 kg) on the comfortnest blue+ model
<b>Flammability</b>	Meets Federal Standard 16 CFR 1632, Cal. TB 117 Optional barrier to meet Cal. TB 129 and Boston Fire IX-11
<b>Guarantee</b>	Non prorated 5-year guarantee on the comfortnest blue and comfortnest blue+ Non prorated 3-year guarantee on the comfortnest blue <sup>MH</sup>

Umano Medical reserves the right to change specifications without notice.



**Umano Medical Inc.**  
230, boul. Nilus-Leclerc  
L'Islet QC G0R 2C0 CANADA  
Toll Free: 1 844 409-4030  
information@umanomedical.com  
umanomedical.com