



UmanoMedical

COMFORT SERIES

comfortnest green™




**Pressure Redistribution • Comfort • Durability**

#### KEY BENEFITS

- 3 layer / 4 zone construction with varying density and IFD ratings for optimal pressure redistribution in the higher risk areas
- Viscoelastic memory foam topper for optimal immersion and comfort
- Entire foam core is constructed of 100% medical grade non-allergenic foam
- All components are anti-microbial and assist in prohibiting growth of bacteria and fungi
- Clam shell design cover features double-folded seams and 4 reinforced handles
- 550 lb weight capacity
- comfortnest green+ model offers High Resilience (HR) foam and a weight capacity of 1000 lb

PRODUCT FEATURES

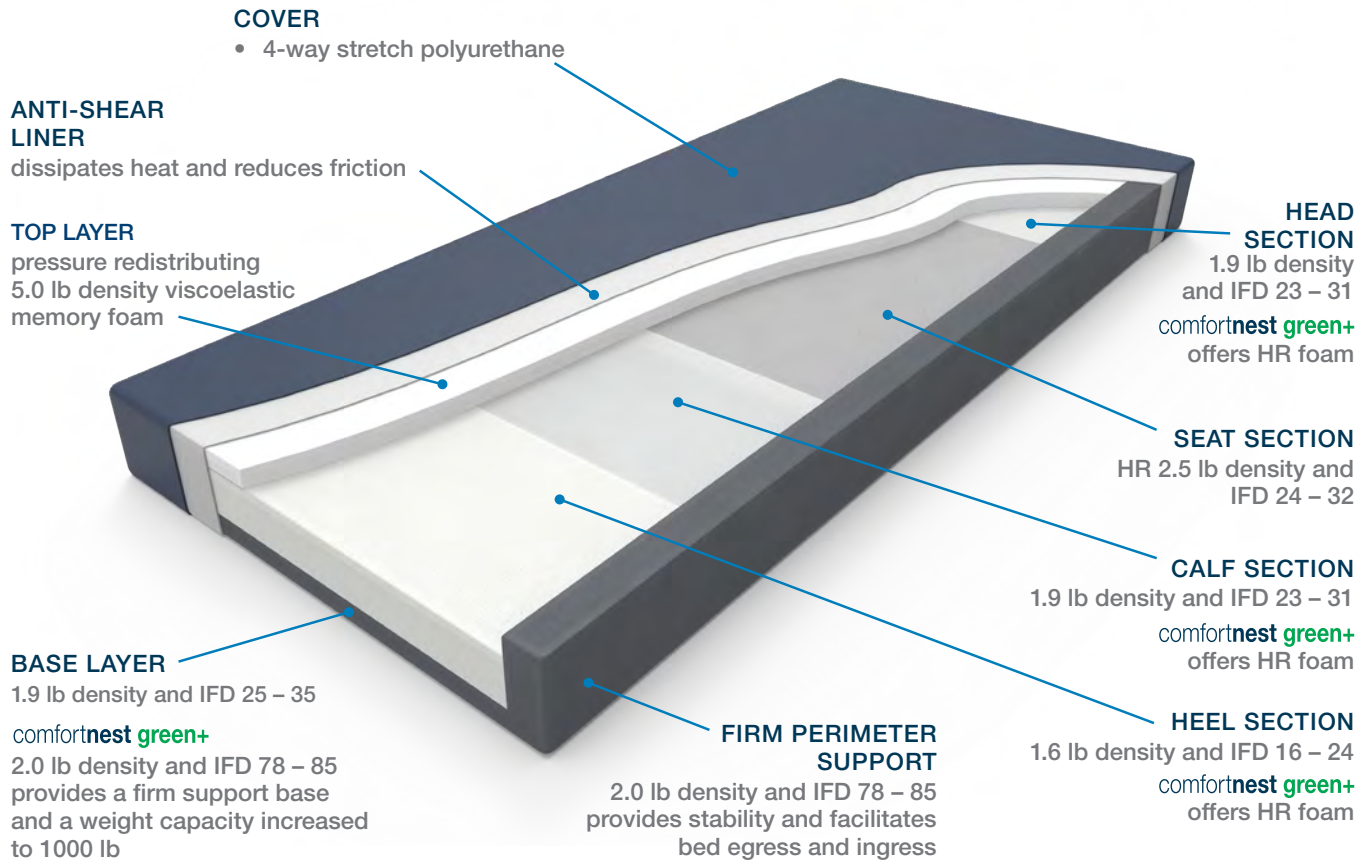

TECHNICAL SPECIFICATIONS

<b>Dimensions</b>	35" W x 80" or 84" L x 6" H (88.9 cm W x 203.2 cm or 213.4 cm L x 15.2 cm H) Custom sizes available
<b>Specialty Model</b>	Behavioral & Mental Health model available: comfortnest green <sup>MH</sup> Fully enclosed cover, no handles, no flaps, no zippers, no anti-shear liner.
<b>Weight Capacity</b>	550 lb (250 Kg) 1000 lb (454 kg) on the comfortnest green+ model
<b>Flammability</b>	Meets Federal Standard 16 CFR 1632, Cal. TB 117 Optional barrier to meet Cal. TB 129 and Boston Fire IX-11
<b>Guarantee</b>	Non prorated 7-year guarantee on the comfortnest green and comfortnest green+ Non prorated 3-year guarantee on the comfortnest green <sup>MH</sup>

Umano Medical reserves the right to change specifications without notice.



**Umano Medical Inc.**  
230, boul. Nilus-Leclerc  
L'Islet QC G0R 2C0 CANADA  
Toll Free: 1 844 409-4030  
information@umanomedical.com  
umanomedical.com