



umano medical

THERAPEUTIC SERIES

covernest gray™



KEY BENEFITS

- Vapor permeable triple layered top cover can be placed over any non-powered pressure redistribution surface to add low air loss therapy.

Therapies

- On-Demand Low Air Loss



PRODUCT FEATURES

CONTROL UNIT

Allows the caregiver to quickly add low air loss therapy to any non powered pressure redistribution surface



TOP LAYER


Is low friction and shear force producing, highly vapor permeable and breathable.

BOTTOM LAYER

is liquid proof to protect the interior of the mattress

MIDDLE LAYER

allows air flow throughout the top cover to assist in prevention of heat build-up and management of skin's microclimate environment.

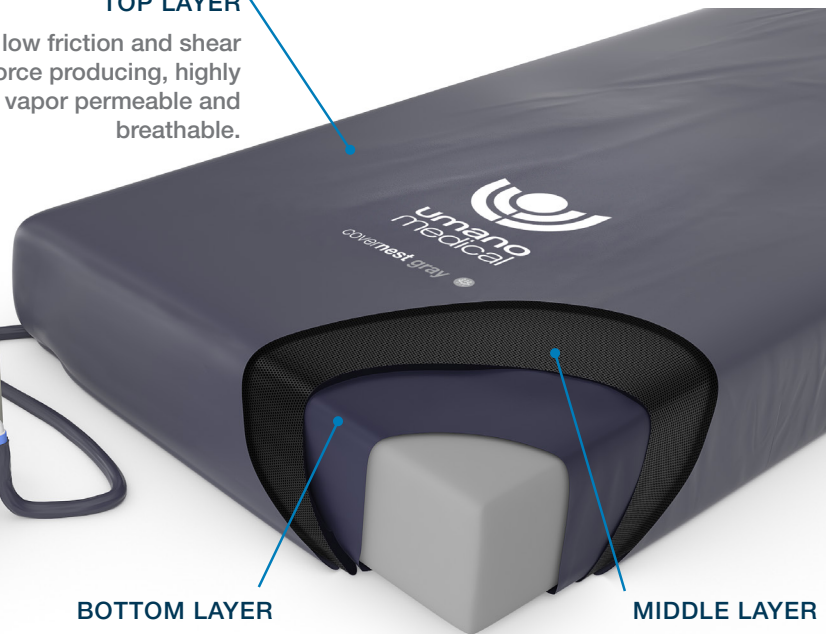

TECHNICAL SPECIFICATIONS

Available Sizes	Width: 36" (91.4 cm)
	Length: 80" and 84" (203 cm and 213 cm)
	Height: Fits 7" to 10" surfaces (17.8 cm to 25.4 cm)
Weight Capacity	as per weight capacity of surface used with LAL Cover
Control Unit	Dimensions 13.8" x 11" x 6.5" (35 cm x 28 cm x 16.4 cm)
	Weight 10.2 lb (4.6 kg)
	Power Cord 14' hospital-grade detachable
	Input Voltage 90 264 V, 50/60 Hz, 1A/0.5A, Universal
	ETL Listed Conforms to UL STD 60601-1 with respect to Electrical Shock, Fire and Mechanical Hazards
Certified to CAN/CSA STD C22.2 No. 601.1	
Flammability	Meets flammability standard CAL 117

Specifications can change without notice.



Manufactured for:
Umano Medical Inc.
 230, boul. Nilus-Leclerc
 L'Islet QC G0R 2C0 CANADA
 Toll Free: 1 (844) 409-4030
 information@umanomedical.com
 umanomedical.com